


## How to Create Dummy CRIS & PAS/HIS Patients via the CRIS System


The CRIS system incorporates utilities for use in Training and Testing including a 'Create Dummy Patient Tool' which enables the creation of Dummy CRIS and PAS/HIS patients.

To access this go to TOOLS > SERVER and select [Create Dummy Patients].

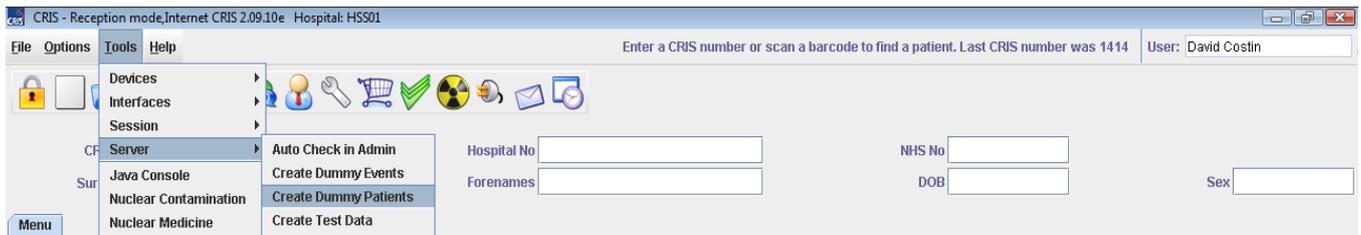

### HOW TO CREATE DUMMY CRIS PATIENTS

Complete the 'Create Dummy Patients' screen as follows:

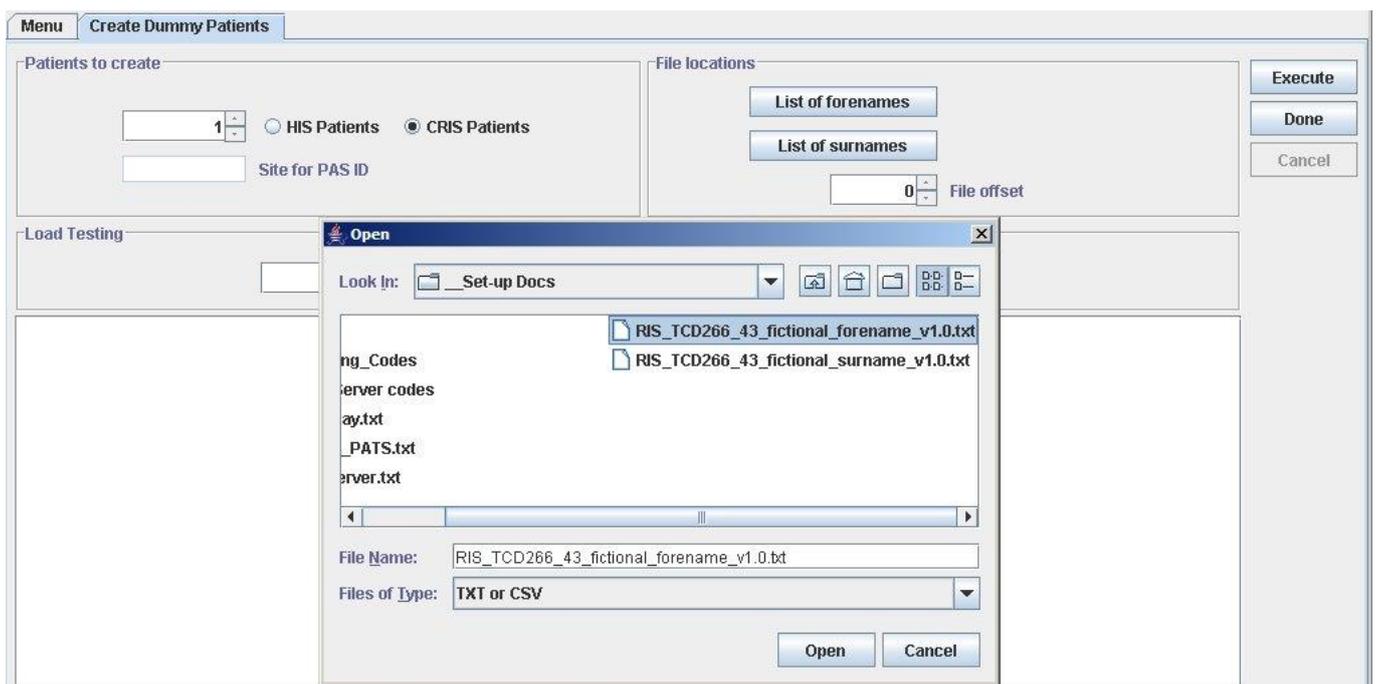

- Click the 'CRIS Patients' Option and enter number of patients contained in the text file (i.e. 43).
- Click [List of forenames] and find and open the appropriate text file – i.e. '43\_fictional\_forename.txt' and click [Open].
- Repeat this process for [List of surnames] and find and opening the text file – i.e. '43\_fictional\_surname.txt' and click [Open].
- Do not amend the file offset.

**Please note:** Increasing the file offset will create additional unique patients but you will need to ensure that you add one extra SURNAME to the 43 Fictional Characters Surname text file. This is because in order for CRIS to randomly generate other fake patients a uneven number of Surnames so that CRIS can step down based on the file offset - i.e. if seeding is 0 it will match Surname One and Forename One but if seeding is 1 it will make Surname Two with Forename One and so on and so forth.

- Click [Execute] to complete the process.
- Once the process is finished you can either proceed to creating ‘Dummy PAS/HIS Patients’ or click [Done].

## HOW TO CREATE DUMMY PAS/HIS PATIENTS

To create Dummy PAS/HIS patients you should repeat the same process as above but this time clicking the ‘HIS Patients’ option and entering the appropriate ‘Site for PAS ID’ before loading the fictional forenames and surname.





- Click [Execute] to complete the process.
- Once the process completes click [Done].

**Please note:** If you are creating dummy patients for use in conjunction with fake orders you need to ensure you link all CRIS patients to PAS patients before proceeding to generate orders.

## FORMAT REQUIRED FOR TEXT FILES

In order to make use of the ‘Create Dummy Patients’ utility it is necessary to create flat text files via notepad using the following format:

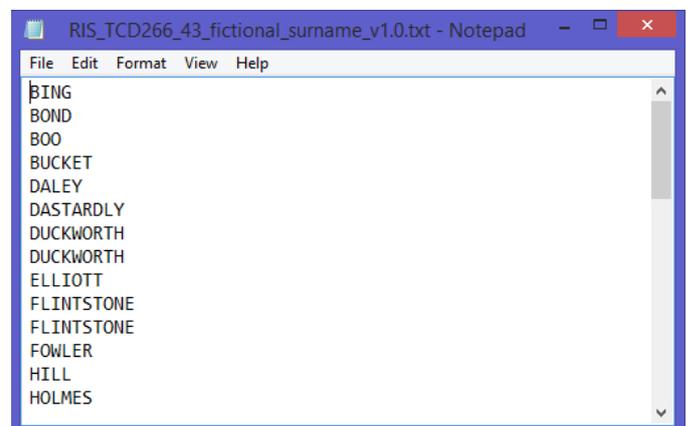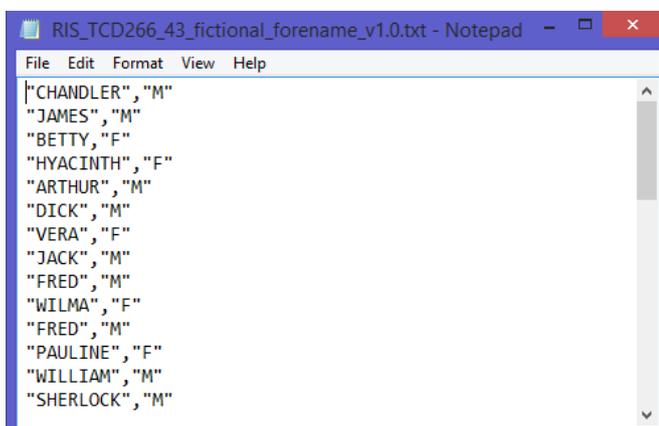

## Security Access Settings

The security setting GENERAL > CREATE\_DUMMY is required via SETUP TABLES > PEOPLE > STAFF > SECURITY SETTINGS in order to access the Create Dummy Patient utility.

## Document Control

<b>Title</b>	How To Create Dummy Patients		
<b>Owner</b>	HSS Training Manager	<b>Date Created</b>	18/06/2014
<b>File Ref.</b>	CRIS_CRIB_CM_286_How_To_Create_Dummy_Patients_V2.0.doc		
<b>CRIS Version</b>	2.09.10p		
<b>Change History</b>			
<b>Issue</b>	<b>Date</b>	<b>Author / Editor</b>	<b>Details of Change</b>
D1.0	18/06/2014	Elena Azzini	Draft Issue
D1.1	08/10/2014	Elena Azzini	Format Changes and New Template
V2.0	29/01/2015	Emma Savage-Mady	Release for 2.09.10p
<b>Review Date</b>	31/01/2016		